# Relativistic Dynamics: The Relations Among Energy, Momentum, and Velocity of Electrons and the Measurement of e/m

MIT Department of Physics

This experiment is a study of the relations between energy, momentum and velocity of relativistic electrons. Using a spherical magnet generating a uniformly vertical magnetic field to accelerate electrons around a circular path and into a narrow slit, you will be able to predict relationships between properties of high velocity electrons. The experiment's goal is to compare your data with the models predicted by nonrelativistic and relativistic dynamics. Additionally you will be able to determine the value of the electron charge to mass ratio, the electron mass, and subsequently the electron charge.

# PREPARATORY QUESTIONS

Please visit the Relativistic Dynamics chapter on the course website to review the background material for this experiment. Answer all questions found in the chapter. Work out the solutions in your laboratory notebook; submit your answers on the course website. [Note: Not available to OCW users]

#### UNITS IN THIS LAB GUIDE

This lab guide uses Gaussian units in which the force on a moving charge in a static field is

$$\vec{F} = q\vec{E} + \frac{q}{c}(\vec{v} \times \vec{B}).$$

The units are:

F — dyne

q — statcoulomb (esu)

B — gauss

E — statvolts cm<sup>-1</sup>

 $v, c - \text{cm sec}^{-1}$ 

Be aware that while magnetometers often display magnetic field values in gauss, voltmeters almost exclusively display voltage in volts, not statvolts.

Some useful constants:

$$e = 4.80298 \times 10^{-10} \text{ statcoulombs}$$
  
 $mc^2 = 511 \text{ keV} = 8.18727 \times 10^{-7} \text{ ergs}$   
1 volt = 3.336 × 10<sup>-3</sup> statvolts

# SUGGESTED PROGRESS CHECK AT END OF $2^{nd}$ SESSION

With the spectrometer magnetic field set to 100 gauss, plot the electron count rate at the detector versus velocity selector voltage (converted to units of  $\beta$ ). Also, plot at least three data points of various kinematic energies versus  $\beta$ .

#### I. INTRODUCTION

Between 1900 and 1910 the relation between the energy, momentum and velocity of charged particles moving at high speeds was a central problem of physics. The fundamental contradictions between Newtonian mechanics and the Maxwell theory of the electromagnetic field, revealed most dramatically in the failures of the Michelson-Morley experiment to detect absolute motion of the earth through the "aether", barred the way to a logically consistent understanding of the deflection of a charged particle by electric and magnetic fields when the particle is moving at a velocity approaching the velocity of electromagnetic waves. Various formulas were derived by Abraham, Lorentz, and Poincare. In 1901 Kaufmann, using the new vacuum techniques pioneered by Thomson, determined the "apparent mass of the electron" by measuring the deflections of the recently discovered  $\beta$  rays from radioactive substances. There was considerable confusion as to whether the experimental results confirmed or contradicted one or another of the formulas. Then Albert Einstein, at the time an obscure 25-year-old examiner in the Swiss patent office, provided a clear and compelling theory of "the electrodynamics of moving bodies" which came to be called the special theory of relativity. Among its remarkable predictions was the slowing of clocks moving at high speed (demonstrated in the Junior Lab experiment on muon decay), and the non-classical relations between momentum, energy and velocity that are demonstrated in the present experiment [1].

In nonrelativistic dynamics a particle acted on by a force  $\vec{F}$  for a time dt over a displacement  $d\vec{r}$  gains momentum  $d\vec{p} = \vec{F}dt$  and kinetic energy  $dK = \vec{F} \cdot d\vec{r}$ .

According to Newtonian dynamics the kinetic energy K, momentum  $\vec{p}$ , and velocity  $\vec{v}$  of a particle are related by the equations

$$\vec{p} = m\vec{v} \tag{1}$$

and

$$K = p^2 / 2m, \tag{2}$$

where  $p^2=\vec{p}\cdot\vec{p}$  and m is the inertial mass of the particle. In nonrelativistic mechanics there is no limit on the

magnitude of  $\vec{v}$ .

According to relativistic dynamics (see [2], or other text on special relativity), these quantities are related by the equations

$$\vec{p} = m\gamma \vec{v} \tag{3}$$

and

$$E^2 = p^2 c^2 + m^2 c^4, (4)$$

where  $E = K + mc^2 = \gamma mc^2$  is the "total relativistic energy." The quantity m is a relativistic invariant and identical with the nonrelativistic inertial mass, and  $\gamma$  is the Lorentz factor defined by the equation

$$\gamma = (1 - \beta^2)^{-1/2},\tag{5}$$

where  $\beta = |\vec{v}|/c$ . Solving equation (4) for the kinetic energy one obtains

$$K = mc^{2} \left( \left[ 1 + \left( \frac{p}{mc} \right)^{2} \right]^{1/2} - 1 \right).$$
 (6)

In the limit of high velocities where  $p \gg mc$ , equation (6) approaches

$$K = pc, (7)$$

which is the exact relation between the energy and momentum of massless particles such as photons and neutrinos. For p < mc equation (6), expanded by the binomial theorem, becomes

$$K = \frac{p^2}{2m} \left[ 1 - \frac{1}{4} \left( \frac{p}{mc} \right)^2 + \dots \right]. \tag{8}$$

In the limit of low velocities ( $|\vec{v}| \ll c$ ) where  $p \ll mc$ , Equation (8) reduces to the nonrelativistic relation expressed by Equation (2).

The electromagnetic force on a charged particle is

$$\vec{F} = q \left[ \vec{E} + \left( \frac{\vec{v}}{c} \right) \times \vec{B} \right], \tag{9}$$

where q is the invariant charge,  $\vec{E}$  and  $\vec{B}$  are the electric and magnetic field intensities, respectively, and c is the invariant speed of light. This force law is valid in both nonrelativistic and relativistic dynamics. In this experiment you will measure the effects of electromagnetic force on the motion of electrons with velocities up to about 0.8c, and measure directly the relations between  $\vec{v}$ ,  $\vec{p}$ , and K.

Some useful references related to beta-decay, beta-spectroscopy and this experiment in particular are given in [3–8].

#### II. EXPERIMENTAL SETUP

The apparatus for the present experiment, shown schematically in Figure 1, is contained in a vacuum chamber inside a spherical magnet that maintains a uniform field in the vertical direction. Inside the vacuum chamber are

- 1. a source of energetic electrons (a minute quantity of  $^{90}$ Sr which emits electrons with a spectrum of energies up to 0.546 MeV yielding a decay product  $^{90}$ Y which emits electrons with a spectrum of energies up to 2.27 MeV),
- 2. baffles that reduce background counts due to scattered electrons that bounce around the vacuum chamber,
- 3. a narrow slit in the baffles at the 90-degree position around a circular path from the source,
- 4. a narrow slit that defines the radius of curvature of the electrons that enter the velocity selector, and
- 5. a solid state PIN diode detector.

The velocity selector has a length of 10 cm and a plate separation distance of  $0.180 \pm 0.003$  cm.

An electron, emitted with a momentum in a narrow range of magnitude and direction, traverses a helical path of fixed radius in the magnetic field and enters the gap between the velocity-selector plates. The distance from the source to the velocity selector, which is also the diameter of the electron's path, is  $40.6\pm0.4$  cm. If the voltage V between the plates is adjusted so that the electric force cancels the magnetic force, then the electron passes through the "velocity selector" in a nearly straight trajectory and strikes the PIN diode. In the diode, a number of silicon valence electrons proportional to the energy deposited by the electron are promoted to the conduction band of the semiconductor and collected as a pulse of charge. The latter is converted by a low-noise preamplifier and precision amplifier into a voltage pulse with an amplitude proportional to the deposited energy. This pulse amplitude is then measured and counted by a multichannel pulse height analyzer (MCA).

For each of several settings of the magnet current one measures

- 1. the magnetic field,
- 2. the voltage across the selector plates that yields the maximum counting rate of pulses in the narrow range of pulse height channels of the MCA corresponding to the energy range of the detected electrons, and
- 3. the median channel of the pulse height distribution.

Given these data, the dimensions of the apparatus, and the energy calibration of the PIN detector, one can determine the momentum, energy, and velocity of electrons and estimate the values of e/m, m, and e.

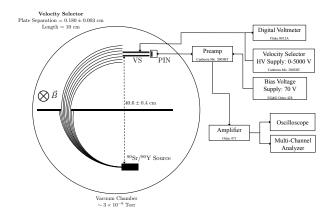


FIG. 1. Schematic diagram of the electron trajectory in the apparatus, the particle spectrometer and associated circuitry. The velocity selector is labeled VS, and the diode detector PIN.

#### II.1. Theory of the Apparatus

In the region between the source and the velocity selector, where only a magnetic field exists, the motion is described by the equation

$$\frac{e|\vec{v}|B}{c} = \left|\frac{d\vec{p}}{dt}\right| = \omega|\vec{p}| = \frac{|\vec{v}|p}{\rho},\tag{10}$$

so

$$B = \left(\frac{c}{e\rho}\right)p,\tag{11}$$

where  $\rho$  is the radius of curvature of the particle trajectory under the influence of the magnetic force. Placement of the source, the collimator, and the aperture of the velocity selector on a circle of radius  $\rho$  allows only particles with a momentum in a narrow range around  $Be\rho/c$  to enter the velocity selector. In the region between the velocity selector plates,  $\vec{E}$ ,  $\vec{B}$ , and  $\vec{v}$  are perpendicular to each other, so one can write the following relation for particles that experience zero deflecting force and go exactly parallel to the plates:

$$eE - \frac{evB}{c} = 0. (12)$$

The voltage between the velocity selector plates is known from Faraday's law,

$$V = \oint \vec{E} \cdot d\vec{\ell},\tag{13}$$

and therefore E = V/d. Hence,

$$\beta = |\vec{v}|/c = E/B. \tag{14}$$

Thus, for any combination of E and B such that E < B, the velocity selector transmits particles with velocities near E/B in a narrow range of magnitudes whose width depends on the geometry of the gap between the plates.

A plot of measured values of B against the ratio E/B reveals the relation between momentum and velocity. According to the nonrelativistic Equation (1), the plot would be fit by a straight line with a slope of  $(mc^2)/(e\rho)$ . Deviation from a straight line as  $E/B \to 1$  indicates the failure of the nonrelativistic relation between momentum and velocity as the velocity approaches c. According to the relativity Equation (3), a plot of B against  $(E/B)[1-(E/B)^2]^{-1/2}$  should be fit by a straight line with a slope of  $(mc^2)/(e\rho)$ . From the slope and knowledge of the values of c and  $\rho$  one can estimate the invariant quantity e/m.

In the experiment it is a good idea to set the magnetic field and then determine the voltage between the selector plates that gives the highest rate of counts of electrons that traverse the circular path and pass between the velocity selector plates to strike the PIN diode detector.

Note that measurements of E and B alone yield a determination of e/m but neither e nor m separately. This is characteristic of all experiments involving only electromagnetic forces. Why is this so? Consider the analogy to the problem of determining the ratio of gravitational to inertial mass of a body moving under the influence of gravity. (Incidentally, the 'Shot Noise' experiment in Junior Lab yields the measurement of e.)

The PIN diode detector combines the virtues of an ultra-thin entrance window and surface dead layer with a total sensitive thickness sufficient to stop electrons with several hundred keV of kinetic energy. Thus, with appropriate calibration, the PIN diode provides a measure of the kinetic energy of the detected electrons. Plots of the kinetic energy against E/B or against  $[1 - (E/B)^2]^{-1/2}$ reveal the relation between energy and velocity; the slope of the latter plot yields a value of m (or more conveniently  $mc^2$  expressed in units of keV in terms of which the energies of the calibration X and gamma-ray photons are expressed). A plot of energy against B reveals the deviation of the energy-momentum relation from the nonrelativistic quadratic form  $E = p^2/(2m)$  toward the linear form E = pc valid for a particle moving with a velocity close to c.

#### II.2. Apparatus Details

The magnet consists of a stack of circular air-core coils enclosing a spherical volume and connected in series so as to produce a current distribution over the surface of the sphere which is approximately equal to the ideal smooth distribution required to produce a uniform field inside the sphere. It turns out that the required distribution of surface current density is

$$J(\theta, \phi) = J_0 \sin \theta, \tag{15}$$

where  $J_0$  is related to the magnitude  $B_i = |\vec{B_i}|$  of the uniform field inside by an equation that is left for the reader to derive (for hints, see Appendix A).

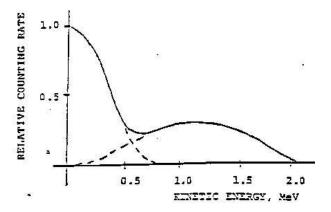


FIG. 2. The  $\beta$ -ray Spectra of  ${}^{90}$ Sr and  ${}^{90}$ Y. The actual energy spectrum of the electrons in this experiment will be distorted by energy loss of the electrons in the material of the source and in the thin metal covering over the source.

The current for the magnet is provided by a stabilized supply whose output can be as high as 200 volts at low impedance, which means that it can deliver a lethal shock. Be careful! Turn off the power supply before touching the magnet. Be sure the cooling fan is turned on whenever the current exceeds 3 amperes.

The coil assembly is split at its equator. The upper hemisphere can be lifted with the block and tackle for access to the spectrometer assembly which is located inside the vacuum vessel under a heavy plastic plate.

The  $\beta$ -ray source,  $^{90}$ Sr, decays with a half life of 28 years to  $^{90}$ Y with emission of an electron with a maximum energy of 0.546 MeV. The yttrium nucleus in turn decays with a half life of 64 hrs to  $^{90}$ Zr with emission of an electron with a maximum energy of 2.27 MeV. In each of these steps, a neutron in the nucleus is transformed into a proton with the simultaneous emission of an electron and an antineutrino. Since the two particles emitted in the  $\beta$ -decay process can share the available energy with the recoiling nucleus in any proportion consistent with the conservation of momentum, the energy spectra of the electrons are continuous from zero to the specified maxima (for details see [4]). A sketch of the energy spectrum of electrons emitted by a thin  $^{90}$ Sr source is displayed in Figure 2.

A narrow slit located at the 90-degree position defines the circular trajectory. The velocity selector consists of two aluminum blocks mounted with their faces precisely parallel to each other and tangent to the central electron path at the entrance to the gap. The voltage for the velocity selector is provided by a high-impedance highvoltage power supply supplemented by a high-voltage battery.

The detector is a PIN diode shielded from light by a thin film of aluminized mylar. The PIN diode should be operated with a positive bias of 50 volts (Caution: voltage above 70 volts will destroy it). A graph of various interaction efficiencies for photons in silicon is

#### Efficiencies in 300 µm Silicon PIN Diode

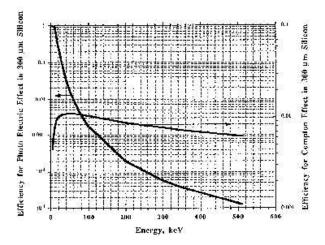


FIG. 3. Photon detection efficiencies in silicon. Note that the vertical axis is logarithmic.

shown in Figure 3.

#### III. EXPERIMENT

Check that the vacuum is better than  $10^{-4}$  torr. A good vacuum is essential because multiple small angle scattering of the electrons by residual gas will distort the trajectories. The rms angle of scattering can be estimated by the following approximate formula:

$$\sqrt{\langle \theta^2 \rangle} = \frac{15 \text{ MeV}/c}{p} \sqrt{\frac{X}{X_0}}, \tag{16}$$

where p is the momentum measured in MeV/c (1 MeV/ $c = 5.34 \times 10^{-17}$  dyne·s), X is the amount of material in the path measured in gm·cm<sup>-2</sup>, and  $X_0$  is a constant, called the "radiation length" and characteristic of the material. The radiation length of air is 37 gm·cm<sup>-2</sup>.

Set up the pulse height measurement chain, and check the performance at each stage with the oscilloscope.

The magnetic field is measured with a Hall effect magnetometer. Before every use of the magnetometer, check the calibration and the zero setting. The Hall sensor, fragile and expensive, is a thin wafer of semiconductor attached to the signal leads (see Figure 4). The wafer and its attached wires should be carefully inserted in the protective aluminum tube and the tube inserted in the plastic block that can be placed on top of the plastic cover of the vacuum chamber. To measure accurately the vertical component of the magnetic field the plane of the semiconductor wafer must be exactly horizontal. To achieve this condition raise the top hemisphere of the magnet just enough to allow you to place the plastic block with the Hall sensor flat inside the aluminum tube in the middle of the vacuum cover where you can assume the field will be almost exactly vertical. Turn on the

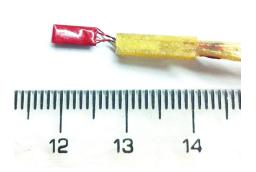


FIG. 4. Hall effect sensor. It is found in at the end of the Hall probe housing.

magnet current to 5 amperes. Rock the plastic block to find the orientation that yields the maximum magnitude of the field. Then rotate the aluminum tube containing the Hall sensor so that the maximum reading is obtained with the plastic block lying flat on the vacuum cover.

Explore the operation of the experiment, and survey the ranges of the measurements before starting a definitive sequence of measurements:

- With the magnet current set to give a field of approximately 100 gauss, measure the precise values of the field at several positions above the circular trajectory from source to PIN diode with the magnet closed to ascertain the degree of uniformity of the field and estimate error in the field measurement due to non-uniformity.
- 2. Calibrate the measurement chain with the X-ray and gamma-ray emissions of the laboratory calibration sources. Be aware of the possibility of a zero offset, *i.e.*, of the 0 channel corresponding to some positive energy. Specially useful for the calibration is the <sup>133</sup>Ba source stored in the lead container under the spherical magnet (<sup>133</sup>Ba decays by electron capture, so in addition to the gamma rays of several energies emitted by the daughter nucleus, the daughter atom, Cs, also emits a K X-ray). You can find the gamma-ray energies by examination of Figure 5 taken from http://ie.lbl.gov/. The K-shell X-ray energies can be found in various X-ray tables (*e.g.* CRC handbook).

The cross section for photoelectric interaction of high-energy photons of energy E in matter of atomic number Z varies approximately as  $Z^4/E^3$ . In the energy range from about 20–100 keV, monoenergetic photons interacting in the silicon PIN diode produce a pulse-height spectrum with a prominent photoelectric peak. Above 200 keV the photoelectric interaction cross section is small relative to the Compton scattering cross section. However, you can use the high end (Compton edge) of

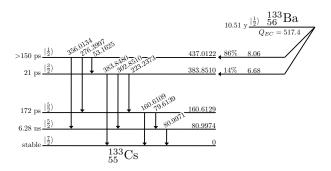


FIG. 5. Decay of  $^{133}$ Ba to  $^{133}$ Cs with associated energy levels. The time on the left is the half-life of the state, next to it is the state's spin, and on the right is the energy in keV. The arrows denote the transitions between energy levels in  $^{133}$ Cs, with energies in keV. Data taken from the Table of Radioactive Isotopes.

the distribution of pulse heights produced by recoil electrons from Compton scattering of the photons as a calibration point. In the text of Melissinos you can find a discussion of Compton scattering, and the relevant information for calculating the energy corresponding to the Compton edge of the distribution of pulse heights produced by electrons that recoil from the scattering events. The first day you should acquire a 2-3 hour (or overnight if possible) calibration in order to generate sufficient events at higher energies for a multipoint calibration. On subsequent days, a 10-15 minute calibration should suffice at the start of each lab session. Plot as many calibration energies as possible against the median MCA channels of the corresponding spectrum features (photopeak or Compton edge). To help correctly identify the relevant peaks and Compton edges consult references [9] and [10]. The resultant calibration should be very linear if you have accurately recognized all the features.

- 3. With the magnet field at approximately 100 gauss, vary the voltage applied to the velocity selector around 4.5 kV while looking for evidence of pulses accumulating in a narrow distribution in the spectrum of pulse heights displayed on the MCA screen. Increase the energy of the selected electrons step by step by increasing B and V, and adjust the measurement chain in particular, the amplifier gain so that the pulses recorded with the highest attainable energies accumulate near the high end of the spectrum displayed on the MCA.
- 4. For each of many (say, 10) values of B spread over as wide a range as possible, measure the voltage V that yields the maximum counting rate. The maximum count rate itself depends on the activity and spectrum of the  $\beta$ -ray source, and on the energy of the incoming electrons, but is typically around tens

of counts per minute. As discussed in the preparatory questions, you can compute a best guess for where in voltage to start searching for the maximum count rate. Scan through a wide enough range of voltages (at least a kilovolt) such that you can establish a statistically significant peak in count rate above background. Make plots in your notebook (by hand, ideally, with error bars) of count rate versus voltage. Do this live in lab as you take the data, so that you may examine your results in real time and adjust your measurement strategy as you go. Keep in mind that the number of counts observed in a fixed time window is a Poisson random variable.

Tabulate B, V, and the kinetic energy corresponding to the median channel number of the pulse-height distribution. Estimate the errors in each of these measurements so that you have the necessary data for error analysis of your results. (Note: the error of a voltage determination is not the FWHM of the measured curve of V vs. count rate.)

## IV. ANALYSIS

For each value of B and  $E = V_0/d$  calculate  $\beta = E/B$ . Then plot  $\beta$  on the y-axis against B on the x-axis. On the same plot draw curves that show the relations between these quantities predicted by nonrelativistic and relativity dynamics.

For each value of B, calculate e/m according to the nonrelativistic and the relativity formulas and plot the results with error bars against B.

Discuss the effect of the spread in the trajectories of the detected electrons on the errors in the determinations of e/m.

Derive a final best value of e/m with a probable error using the techniques in [11].

Plot your measured values of K against p, and draw the curves predicted by nonrelativistic and relativity mechanics.

#### V. POSSIBLE THEORETICAL TOPICS

- 1. The surface current distribution required to produce a uniform field inside the spherical magnet.
- 2. Relativistic dynamics *i.e.*, the relation between v, p, and E.
- 3. Compton scattering and the photoelectric effect at gamma-ray energies.
- 4. PIN diodes as particle detectors.
- 5. Resolution of the velocity selector.
- 6. Multiple Coulomb scattering of the electrons by the residual air in the vacuum chamber.
- [1] Various, selected reprints by the American Institute of Physics (1963).
- [2] A. P. French, Special Relativity (MIT Press, 1968).
- [3] E. Bleuler and G. Goldsmith, *Experimental Nucleonics* (Holt, Rinehart, and Winston, 1960) pp. 342–346.
- [4] R. Evans, The Atomic Nucleus (McGraw-Hill, 1953) chapter 17, pp. 536–566, Also contains statistics in chapter 26.
- [5] J. England, Techniques in Nuclear Structure Physics, Vol. Vol. II (Wiley and Sons, 1974) pp. 349 ff.
- [6] R. B. L. R. P. Feynman and M. Sands, The Feynman Lectures on Physics (Addison-Wesley, 1963) i-15.
- [7] J. D. Jackson, *Classical Electrodynamics* (Wiley and Sons, 1962).
- [8] D. J. Griffiths, Introduction to Electrodynamics (Upper Saddler River, 1999).
- [9] M. McDonnell and M. K. Ramaswamy, Phys. Rev. 168, 1393 (Apr 1968).
- [10] B. Crasemann, J. G. Pengra, and I. E. Lindstrom, Phys. Rev. 108, 1500 (Dec 1957).
- [11] P. R. Bevington and D. K. Robinson, Data Reduction and Error Analysis for the Physical Sciences, 3rd ed. (McGraw-Hill, 2003).

#### Appendix A: The Spherical Uniform-Field Magnet

The magnet is an arrangement of currents in coaxial coils distributed over the surface of a sphere in such a way as to produce a uniform magnetic field  $B_i$  inside the sphere. The design problem is to find the distribution of surface current density  $J(\theta,\phi)$  that achieves this condition.

Everywhere but on the surface of the sphere  $\nabla \times B = j = 0$ . It follows that  $B = -\nabla \phi_m$ , where  $\phi_m$  is the magnetic potential function, *i.e.* a solution of Laplace's equation  $\nabla \times (\nabla \phi_m) = \nabla^2 \phi_m = 0$ . Potential theory says that a solution of Laplace's equation that satisfies a particular set of boundary conditions, such as the distribution of currents on the surface of our sphere, is unique except for an additive constant. Thus if you can cook up a solution of Laplace's equation that does, indeed, satisfy the boundary conditions, you can be sure you have got the one and only solution to the problem.

If the field is to be uniform inside the sphere, then

$$\phi_m(r, \theta, \phi) = -B_i z = -B_i r \cos \theta.$$
  $r < R$  (A1)

Outside at large distances it is certainly true that

$$\phi_m = \frac{m\cos\theta}{r^2} \qquad r \gg R,$$
 (A2)

i.e. the potential must have the form of a dipole potential at large distances since the field is the sum of contributions from many coaxial circular current loops. The problem boils down to that of matching an outside dipole field to the inside uniform field  $B_i$  by proper adjustment of the surface current distribution in such a way as to satisfy the continuity requirements on the radial and tangential components of  $\vec{B}$ . The continuity conditions follow directly from the two Maxwell equations governing the magnetostatic field, namely

$$\vec{\nabla} \cdot \vec{B} = 0,$$
 
$$\vec{\nabla} \times \vec{B} = \left(\frac{4\pi}{c}\right) \vec{j}.$$

(Note:  $\vec{j}$  is the volume current density with dimensions of statamp·cm<sup>-2</sup>, whereas what we seek is the surface current density with dimensions statamp·cm<sup>-1</sup>.)

The rest of the solution is left to the reader.

For the general approach to the solution of such boundary condition problems see, e.g., Jackson [7] or Griffiths [8].

## Appendix B: Equipment List

Manufacturer	Description	URL
Agilent	Oscilloscope	agilent.com
Sorenson	0–5A Power Supply	
RFL	Gaussmeter	
Bertran	7.5kV Velocity Selector PS	
Canberra 2006	Charged Particle Preamp	
Hamamatsu	PIN Diode	
Ortec	Amplifier	ortec-online.com

MIT OpenCourseWare https://ocw.mit.edu

8.13-14 Experimental Physics I & II "Junior Lab" Fall 2016 - Spring 2017

For information about citing these materials or our Terms of Use, visit: <a href="https://ocw.mit.edu/terms">https://ocw.mit.edu/terms</a>.